



*Picture of the model with optionnal headrest cushion ACCCOUSSIN-01  
option 4 wheels with central brake system ACC150-MULT and  
ACCMUL-RPIEDS Electric foot rest with cushion*

*Electric multi-treatment chair*

**MULTIKA version V2 range**

**MULT-V2-65-ME2**

**CARINA**  
Medical




**CE**  
ISO 9001  
ISO 13485

## TECHNICAL DETAILS

- 3 sections multi-treatment chair, width 65cm
- Control through multifunction hand remote with buttons for memory positions ( Seat position and transfer position )
- Adjustable backrest and legrest with electric jacks
- Tredelenburg position
- One section mattress , thickness 90mm (6cm of HR35 foam + 3cm of viscoelastic memory foam)
- Equipped with a pair of armrest on knee system, synchronized with the backrest, a push bar into the backrest, and a pair of DIN rails 25x10 to fix accessories
- Shapped base for ideal cleaning, base colour light grey (RAL 7035)
- On 4 design wheels with one brake on each wheel diam 125mm

## DIMENSIONS

- Width 65cm, width overall 85cm
- Adjustable height from 56 to 96cm
- Load capacity: 200kg in dynamic

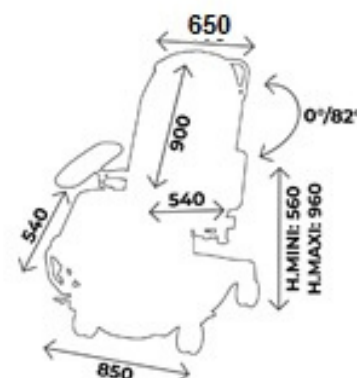

### You may also like :

Réf. **MULT-V2-75-ME2**: Same as MULT-V2-65-ME2 but width 75cm



### WHAT YOU NEED TO KNOW

- ✓ Delivery time: 15 days upon receipt of order
- ✓ Warranty: 5 years (1 year for the upholstery)

# ZOOM MULT-V2-65-ME2 with accessories




# ACCESSORIES



## ACC MUL-BPPIEDS

Retractable push handles on the foot side

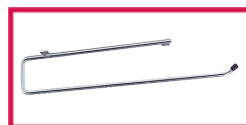

## ACC 125

Paper roll holder for exam sheet



## ACC BATTERIE01

Battery, capacity 20 cycles  
(1 cycle = 1x height's up/ down +  
1x backrest's up/down)



## ACC COUSSIN-01

Comfortable cushion  
50 x 40cm



## ACC CLAMEDIAM16

Adaptation clamp in aluminium  
for fixation on rail dim 25x10, the unit  
siam.14mm, diam.16mm



## ACC MULT-320INO1

Adjustable infusion stand with 2 plastic  
hooks, without clamp



## ACC 514-MULT65ME

Plastic protection for leg rest



## ACC MULTIK08-AM

Oxygen bottle holder with fixation  
straps



## ACC 150-MULT

4 wheels with central brake system  
diam 150mm



## ACC-MUL-RPIEDS

Electric foot rest with cushion



## ACC MULTIK02-PIED

Emergency pedal to be put in  
Trendelenburg position quickly



## ACC EQUIPOTENT01

Equipotentiality plug



## ACC MUL-5ROUES

5th wheel to optimize transfers  
and mobility of the chair, to be  
order with ACC150-MULT option



## ACC ROUE-MOTOR01

Motorized assistance



## ACC MULTIK22


Manual release of the backrest

# UPHOLSTERY'S COLOURS

Fire retardant  
upholstery, M1 class

*For the upholstery colours, you can choose between the Classic and Urban range*

## *Upholstery : Classic range*

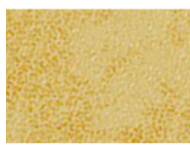

**Arno**



**Garonne**




**Potomac**




**Amour**




**Mississippi**

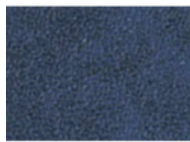


**St Laurent**

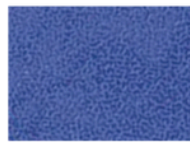

**Rhône**

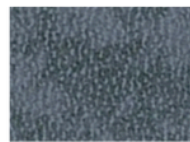

**Missouri**

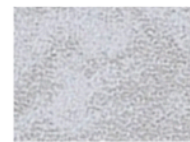

**Seine**



**Amazone**

## *Upholstery : Urban range*



**Abysses**



**Cendre**



**Origan**




**Géranium**




**Réglice**




**Prune**

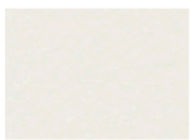


**Perle**

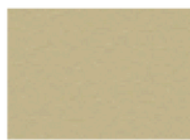


**Ficelle**

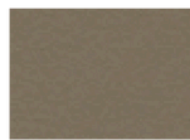


**Terre**

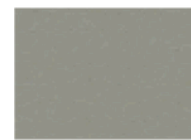


**Béton**

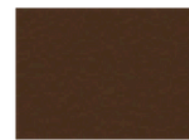

**Chataigne**

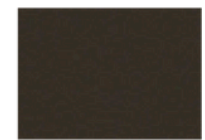

**Chocolat**

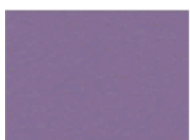

**Pensée**




**Pétrole**



**Turquoise**



**Granny**

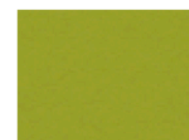

**Absinthe**



**Ananas**



**Orange**



**Fuschia**



**Tomate**



# MAINTENANCE GUIDE

## **DAILY MAINTENANCE**

Upholstery : Cleaning with soapy water only (no solvents), can be disinfected with products such as Phagosept Spray, Salvanios pH10 at 0.75% and Anios Surfaces Hautes, without abrasive material

Base : cleaning with soapy water only (no solvents) without abrasive materials

## **WARRANTY**

Mechanical/ electrical part: 5 years

Upholstery: 1 year



CONTACT

CARINA  
*Médical*

Parc d'entreprises Visionis  
278 rue de l'industrie  
01090 GUÉREINS- FRANCE

**COMMERCIAL SERVICE**

Tel: +33(0)4 74 66 69 49  
Fax: +33(0)4 74 66 69 50  
[export@carinamedical.com](mailto:export@carinamedical.com)