



Picture of the model with options ACCCOUSSIN-01: Comfortable cushion, and ACCMUL-BPPIEDS retractable pushing handles on foot side

Electric multi-treatment chair

MULTIKA version V2 range

MULT-V2-75-AM2

CARINA
Medical




CE
ISO 9001
ISO 13485

TECHNICAL DETAILS

- 3 sections multi-treatment chair, width 75cm
- Control through multifunction hand remote control with buttons for memory positions (Seat position and transfer position)
- Adjustable backrest and legrest with electric jacks
- Tredelenburg position
- One section mattress , thickness 90mm (60mm of HR35 foam + 30mm of viscoelastic memory foam)
- Pair of retractable security side rails with ergonomic and safe handles, and a fix pushing bar on the head side. The side rails are designed to be able to support an optionnal infusion stand
- Shapped base for ideal cleaning, base colour light grey (RAL 7035)
- With 4 wheels diam.150mm with central brake system: Two pedals at the front of the chair and a brake bar on the head side

DIMENSIONS

- Width 75cm, width overall 85cm
- Adjustable height from 56 to 96cm
- Load capacity: 200kg in dynamic

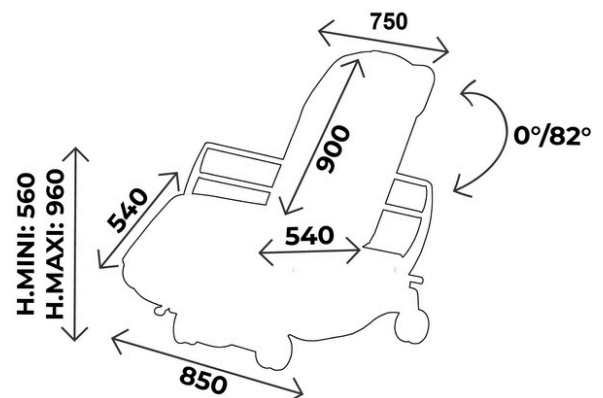

You may also like :

Réf. **MULT-V2-65-AM2**: Same as MULT-V2-75-AM2 but width 65cm



WHAT YOU NEED TO KNOW

- ✓ Delivery time: 15 days upon receipt of order
- ✓ Warranty: 5 years (1 year for the upholstery)

ZOOM
MULT-V2-75-AM2
With accessories

10 buttons multi-function hand remote
Fixation on the left or on the right side

ACCCOUSSIN-01 (option)

Comfortable cushion 50 x 40cm

One section mattress

Upholstery with visco-elastic foam for an optimal comfort

Pair of security side rails with support to fix infusion stand

ACCMUL-BPPIEDS (option)

Retractable and adjustable pushing handles

Push handle

Electric backrest and leg rest


Braking bar on the back side

4 wheels diam 150mm with central brake system, two pedals in front of the couch

ACCESSORIES



ACC MUL-BPPIEDS
Retractable push handles on the foot side

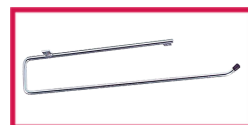

ACC 125
Paper roll holder for exam sheet



ACC MULT-320INO1
Adjustable infusion stand with 2 plastic hooks



ACC COUSSIN-01
Comfortable cushion
50 x 40cm



ACC 514-MULT75AM
Plastic protection for leg rest



ACC MULTIK05-BA2
Reading tray



ACC MULTIK30
Holder for medical files, to be fixed into the backrest



ACC MULTIK08-AM
Oxygen bottle holder with fixation straps



ACC MULTIK02-PIED
Emergency pedal to be put in Trendelenburg position quickly



ACC BATTERIE01
Rechargeable battery to power the chair when not connected to main socket



ACC EQUIPOTENT01
Equipotentiality plug



ACC MUL-5ROUES
5th wheel to optimize transfers and mobility of the chair



ACC ROUE-MOTOR01
Motorized assistance




ACC MUL-RPIEDS
Electric foot rest, (not compatible with ACCMUL-BPPIEDS: pushing handles on foot side)

UPHOLSTERY'S COLOURS

Fire retardant
upholstery, M1 class

*For the upholstery colours, you can choose between the Classic and Urban range
Be careful ! Urban colours with extra value*

Upholstery : Classic range

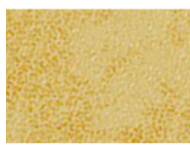

Arno



Garonne




Potomac




Amour




Mississippi

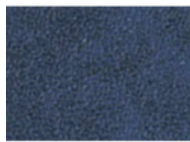


St Laurent

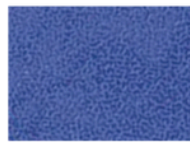

Rhône

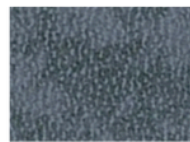

Missouri

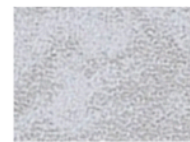

Seine



Amazone

Upholstery : Urban range



Abysses



Cendre



Origan




Géranium




Régilisse




Prune

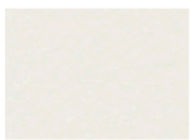


Perle

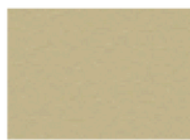


Ficelle

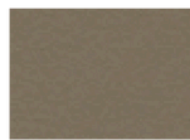


Terre

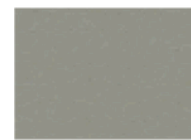


Béton

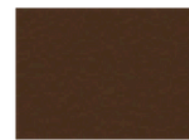

Chataigne

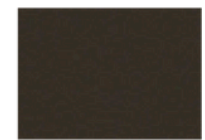

Chocolat

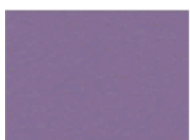

Pensée




Pétrole



Turquoise



Granny

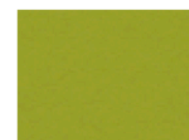

Absinthe



Ananas



Orange



Fuschia



Tomate



**BASE
COLOURS**



RAL 9003



RAL 1013



RAL 1015



RAL 6019



RAL 6033



RAL 6000




RAL 5017



RAL 5002



RAL 5000

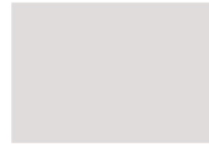

RAL 7035



RAL 7016



RAL 9005




RAL 1018



RAL 3020



RAL 3003

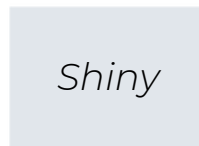

SILVER



MAINTENANCE GUIDE

DAILY MAINTENANCE

Upholstery : Cleaning with soapy water only (no solvents), can be disinfected with products such as Phagosept Spray, Salvaniol pH10 at 0.75% and Anios Surfaces Hautes, without abrasive material

Base : cleaning with soapy water only (no solvents) without abrasive materials

WARRANTY

Mechanical/ electrical part: 5 years

Upholstery: 1 year



CONTACT

CARINA
Médical

Parc d'entreprises Visionis
278 rue de l'industrie
01090 GUÉREINS- FRANCE

COMMERCIAL SERVICE

Tel: +33(0)4 74 66 69 49
Fax: +33(0)4 74 66 69 50
export@carinamedical.com